
Differentiating **iPSCs**

Which approach works best?

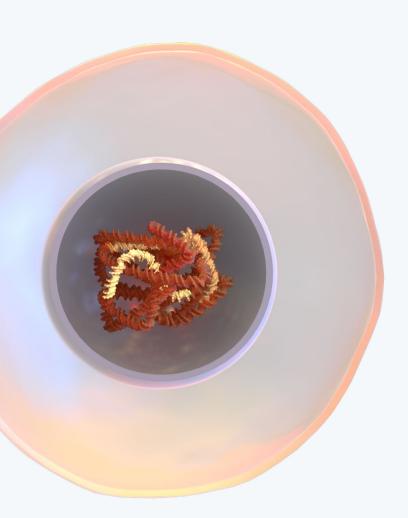

Overcoming the challenges of cell differentiation and reprogramming with opti-ox™, precision reprogramming

Current methods to generate cells from human induced Pluripotent Stem Cells (iPSCs) - directed differentiation and cellular reprogramming - lack reproducibility and scalability. This raises challenges for disease research, drug discovery and cell therapy development due to unreliable data and long experimental timelines.


A technology that enables the reproducible, consistent and scalable production of physiologically relevant human cell models for any cell type would mean scientists can be more confident in their experiments. opti-ox, a first-of-its-kind precision reprogramming technology, meets all of these requirements by controlling the precise expression of cell-fate determining transcription factors without gene silencing.

In this infographic we outline traditional methods to generate cells from human iPSCs and compare these to the new precision reprogramming technology enabled by opti-ox.

Directed differentiation


- Low purity
- Low scalability
- Inconsistent batches
- Long laborious complex protocols

Directed differentiation generates somatic cells from iPSCs through extracellular signalling by the addition of growth factors, cytokines or small molecules mimicking developmental signals, driving cell fate.

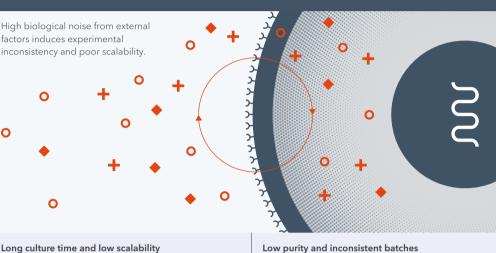

Long, complicated protocols combined

with the costly use of reagents (growth factors, proteins etc.) limit the ability to scale.

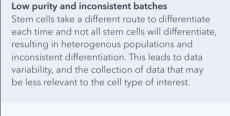






Cellular reprogramming

involves the ectopic expression of transcription factors (TFs) to drive a specific cell fate. With this method cell identity is less dependent on external cues and has improved cell specificity. However the random integration of TFs can result in gene silencing meaning cell populations are inconsistent between batches.

Cellular reprogramming

TFs drives a defined cell fate but is susceptible to gene silencing.

Reduced biological noise from

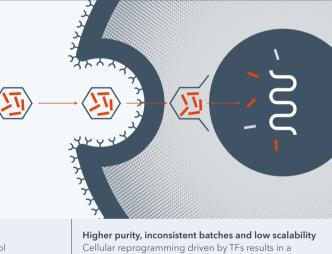

more homogenous population compared to directed

differentiation. However the random integration of TFs leads to variable silencing resulting in inconsistent reprogramming, heterogeneity across batches and

Low scalability

Fast protocol

Inconsistent batches

Cell reprogramming offers a faster protocol with differentiation occurring over weeks rather than months.

Shorter culture time



opti-ox[™] precision reprogramming by bit.bio

Consistent batches

Highly pure with defined identity


poor scalability.

a controlled gene engineering approach that enables unlimited batches of any human cell to be manufactured consistently at scale. opti-ox utilises two genome safe harbour sites to protect TFs

reprogramming technology is

opti-ox™ precision

from silencing, and leverages an inducible system to control TF activation. With TFs precisely switched on in every stem cell, entire cultures consistently and rapidly convert into a defined cell type on a vast scale.

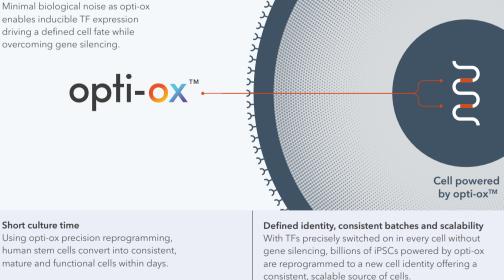



bit.bio combines the concepts of cell programming and biology

at scale.

Who we are

to provide human cells for research, drug discovery and cell therapy, enabling a new generation of medicines. This is possible with our precision human cellular reprogramming

technology opti-ox™ - a gene engineering approach that enables

unlimited batches of any human cell to be manufactured consistently

To order or speak with a member of our team, email orders@bit.bio.

bit.bio

Contact us

The Dorothy Hodgkin Building Babraham Research Campus Cambridge CB22 3FH United Kingdom +44 (0) 1223 787 297

www.bit.bio

15% off your first vial. Are you interested in testing how bit.bio's precision reprogrammed ioCells could advance your experiments?



Use code **OPTIOX15** when speaking to your account manager for 15% off your first vial. Terms and conditions apply*.

* Terms and Conditions: Use of the promotion code is required to receive your discount. This offer is one time use only and cannot be combined with any other price discounts, volume purchase agreements, university programs, standing order discounts or other contract pricing. Provide this code to your bit bio representative when you place your next order. Discount valid only on Small and Large vials of ioGlutamatergic Neurons (io1001), ioSkeletal Myocytes (io1002) and ioGlutamatergic Neurons HTT 50CAG/WT (EAio1004). Discount cannot be applied to prior orders or purchases. Other restrictions may apply. Void where prohibited. Offer expires 31 October 2022.



THE CELL CODING COMPANY